

STEEMA

DUAL TOT SLING SYSTEM

Transobturator sling systems comparison


	TVT-O [Gynecare]	Obtryx II [Boston Scientific]	STEEMA [Promedon]
Mesh Thickness	0.63 mm ¹	0.66 mm ¹	0.5 mm ³
Fiber size [diameter]	0.15 mm ¹	0.15 mm ¹	0.14 mm ³
Weight [g/m ²]	100 ¹	100 ¹	90 ³
Load at failure	70 N ¹	70 N ¹	90 N ³
Sling delivery force	1,11 Kg ⁴	0,73 Kg ²	0,77 Kg ⁴
Tunneling force	1,65 ⁴ Nm	1,16 ⁴ Nm [Halo needle]	1,15 ⁴ Nm

ORDERING INFORMATION

STEEMA Dual TOT Sling System

Product code: KIT-STEEMA-T PLUS

- 1 | STEEMA Sling
- 2 | Needles STEEMA
- 1 | Urethral Protector

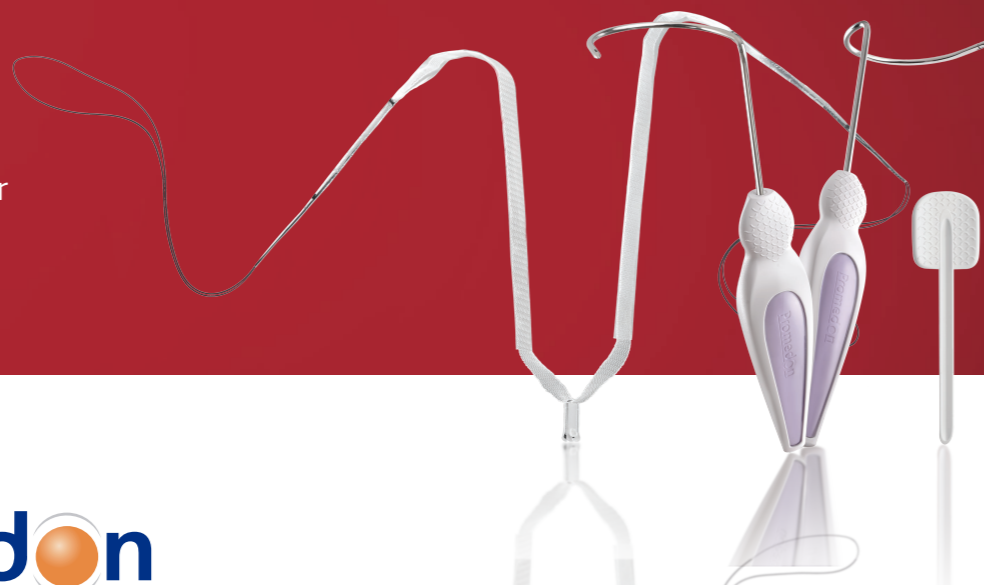

CE 0197 B-30-76(02)/10-FEB-2018


Promedon
People + Innovation

www.promedonproducts.com/steema

1. Moali, P., et al. Tensile properties of five commonly used mid-urethral slings relative to the TVT™ - Int Urogynecol J [2008] 19:655-663.
 2. Information published by Boston Scientific: Ref.: WH-119702-AA XM 11/12.
 3. Data on file. Promedon sling design properties.
 4. Data on file. Bench test sample size n=20. Results from these studies may not be predictive of results in other cases.

STEEMA

DUAL TOT SLING SYSTEM

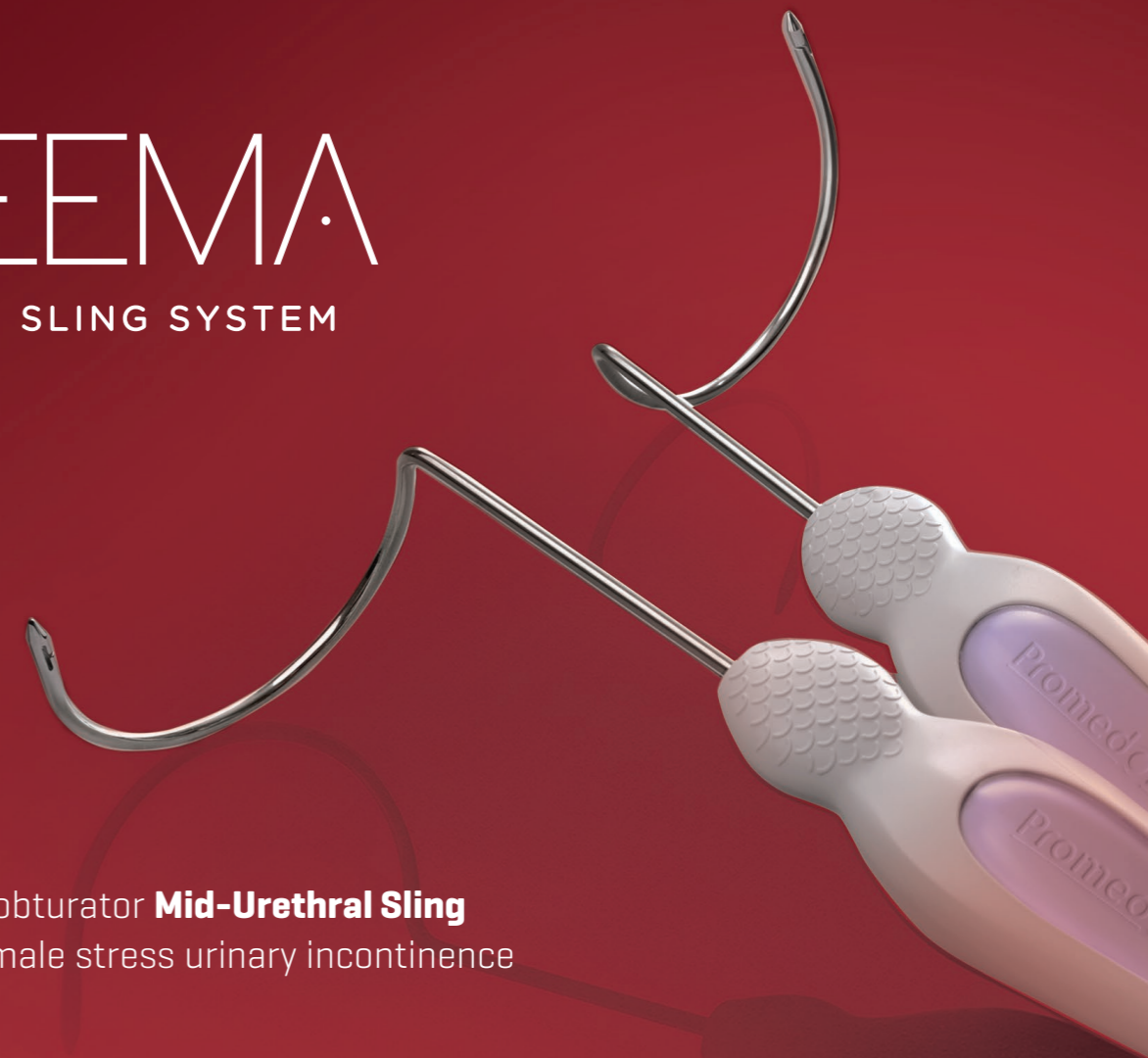

Premium Transobturator **Mid-Urethral Sling** system to treat female stress urinary incontinence



DUAL APPROACH

EASE OF POSITIONING

LOW ELASTICITY

Promedon
People + Innovation

BRINGING THE ULTIMATE VERSATILITY IN TOT PROCEDURES

Premium Transobturator Sling System
for SUI treatment. Designed to be used
either in-out or out-in, depending on
surgeon's preference.



Reduced Needle Diameter

Smallest needle diameter (3mm)
on the market causing less tissue
disruption and allowing a smooth
passage.

Handle


Ergonomic design ensures
precise, controlled maneuvering.

T-Shaped
eyelet

HIGHLIGHT

**30% LESS
TUNNELING FORCE**

COMPARED TO THE CATEGORY LEADER ¹

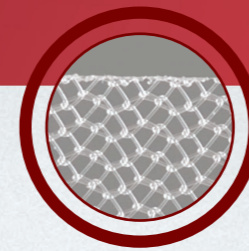

Polypropylene type I mesh

Thermosealed edges provide
better structural support.
Low elasticity mesh.



Dilator

Ensures a smoother
sling passage.




Protective Sheath

Minimizes friction and tissue
damage when passing
the sling.



Positioning Tab

Provides a controlled
and easy positioning of
the sling. Removable.

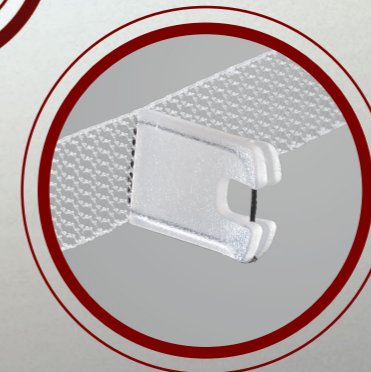

Urethral Protector

Ergonomic grip.
Insertion mark.



1. DATA ON FILE. BENCH TEST SAMPLE SIZE N=10. RESULTS FROM THESE STUDIES MAY NOT BE PREDICTIVE OF RESULTS IN OTHER CASES. PRODUCTS TESTED STEEMA VS THE IN-OUT MARKET LEADER.